

**VEGAS**  
COSMETICS

**SKIN CELL TRONIC 2.0**

› User Manual



› Photon ultrasonic anti aging device

## Instructions for Use of the Device

- › Thank you very much for choosing our Skin Cell Tronic device.
- › Please read the instructions carefully, so you will be adequately informed before you start using this device. Please keep the instructions in a safe place for future reference.
- › This device was designed for the sole purpose of facilitating beauty treatments in the privacy of your own home. Use of this device to treat any medical condition or illness of any kind is not permitted and the manufacturer will not be liable for any consequences arising out of utilization methods other than explicitly suggested in the user manual.

### Contents

◆ Important information regarding Skin Cell Tronic	1
◆ Safety information	2
◆ Functions information	4
◆ Charging your device	6
◆ Skin care instructions	7
◆ General care instructions	7
◆ Anti aging care concept	8
◆ Even skin tone care concept	10
◆ Acne and sensitive skin care concept	10
◆ Warranty card	12

### Technical Data

◆ Product name:	Skin Cell Tronic
◆ Model:	SCT02 – Product code: 460
◆ Adaptor input power:	AC 100-240V
◆ Adaptor output power:	15V, 400 mA
◆ Dimensions:	170 x 46 x 50mm (6.7 x 1.8 x 2 in)
◆ Weight:	132g (4.7 oz)
◆ Ultrasound frequency:	3MHz±5%



## Important Information regarding Skin Cell Tronic

- ◆ This ultrasound skin care device delivers up to three (3) million vibrations per second. It combines a high level of modern technology and efficiency with utmost safety and peace of mind. As the reach of the ultrasound vibrations is limited to the upper layers of the skin, its use is most suited for the eye area. High ultrasound frequencies stimulate the skin cells in a comfortable manner and help activate the cell metabolism.
- ◆ The light treatment incorporates the use of different colored lights, each with specific functions in the treatment process:
- ◆ The blue light acts in a calming manner and is recommended for problem skin (i.e. rashes, sensitive skin, oily skin with blemishes). It can help adjust oil production to a more normal level, effectively eliminating acne and blemishes, as well as tighten the skin in order to achieve a healthy and beautiful appearance.
- ◆ The red light is suitable for all skin types. Skin cells are being stimulated by the warmth of the light treatment which results in improved blood circulation and increased production of collagen and protein. As your skin becomes more adaptable you will notice a tighter and more youthful appearance. The mildness of the red light wave length is especially appropriate for use on forehead and around the sensitive eye area. The green light is most suitable for combination skin and activates micro-circulation as well as increases oxygen intake with calming and relaxing effects.
- ◆ The different colored light waves accelerate moisture absorption from the skin care products used during and immediately after the treatment process.
- ◆ The overlapping effects of different colored lights (blue, red and green) combined with the 3M ultrasound waves and the ionic function make this skin care device truly remarkable.
- ◆ Our ultrasound device was carefully designed to adapt well to the facial contours. You will find it easy to use and operate.

## Safety Information

- › This product is a device specifically designed to beautify healthy skin, and explicitly disallowed for any other use.
- › Carefully observe the manufacturer's instructions regarding frequency of use and functions. We strongly suggest limiting the very first treatment time exposure. This is especially important if your skin tends towards allergic reactions.
- › Immediately stop using this device if you experience any discomfort.

## Under no circumstances should you use this device if

- › You are pregnant
- › You have suffered facial fractures and/or had facial reconstructive surgery or any kind of metal, silicone or plastic facial implants
- › You have sunburn or skin irritation due to excessive sun exposure
- › You experience skin irritation caused by cosmetics
- › You suffer from dermatitis
- › You experience itching and elevated skin temperature
- › Your skin is irritated due to a physical condition
- › You are suffering from enlarged arterioles (small blood vessels) in your face caused by long-term use of steroids, hormones or due to liver dysfunction
- › You have metal implants anywhere in your body
- › You suffer from heart disease (particularly patients with a pacemaker/defibrillator)
- › You suffer from any condition not specifically mentioned in this operation manual

**Never use this device directly on your throat or on your eyes!**

## Safety Instructions

In the interest of your own personal safety, we urge you to read the following information carefully before using this cosmetic device. Our product is a well-designed precision instrument. In order to achieve perfect and satisfactory performance, please read and follow the instructions outlined below with care. This will not only ensure a safe and successful treatment process but also secure the full extent and validity of the product warranty.



## An easy test for potential allergies

1. Cleanse your wrist thoroughly.
2. Apply cream (we recommend Hyaluronan Gel), leave on your wrist and do not rinse off for 48 hours. Stop test immediately should you experience any discomfort during this test period.
3. If you do not experience any negative reactions during the test period, you may proceed to use the device.
4. Stop using this product immediately and consult a physician if you experience any sign of discomfort or an allergic reaction.



## Operations and Maintenance Precautions

1. Store product in a safe, dry place, do not use under water. Otherwise, the electronic components could become corroded or short circuit.
2. Never attempt to open or disassemble the product yourself. All product components must be maintained and/or repaired by a trained and qualified professional technician. Mishandling the device will result in the warranty being voided!
3. Do not use any adaptor other than the one included with this product in the original packaging. Using an incompatible adaptor may result in malfunction and will also void the warranty.
4. Keep this product away from children as they may accidentally swallow small parts.



WARNING!

### Important Safety Notice



Disassembling the product yourself is not permitted



Keep this product away from children




Do not expose to extreme heat or fire



Store product in dry places and do not use under water



Do not drop or set down roughly

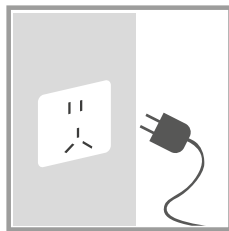

Always unplug device after use

## Functions Information

- Ultrasound with up to three (3) million vibrations per second
- Ionic function (positive and negative)
- Contact surface made out of high-quality stainless steel

Red LED light: 12 LED  
Blue LED light: 12 LED  
Green LED light: 12 LED  
High energy, fully focused



Ionic electrode bar



Message intensity indicator

Ionic+ or Ionic- working indicator

Ultrasound working indicator

Power/charging indicator

Pre-set treatment programmes

Ionic+ or Ionic- selection

Sonic (ultrasound) selection on/off

Photon selection (LED light)

Power on/off

CLEANS

ANTI-AGE

SENSITIVE

ION

SONIC

POWER  
CHARGE

COMBI

ION

SONIC

PHOTON

ON/OFF




- You need to touch the ionic electrode bar with your finger while using the ionic function in order to keep the power circuit closed. If you are using the device on somebody else, you have to stay in physical touch with that person.

## Charging your device

- The “power on” light will blink as long as the device is charging. The light will go off as soon as the charging process is complete.
- As soon as the charge drops below 3V, all lights will blink 10 times and the device will shut off automatically.
- Your Skin Cell Tronic device comes with a rechargeable lithium battery. Please make sure that you fully recharge it every time it runs down (see above).
- You may use the device even while you are charging it. Before using it for the first time, the device should be fully charged.

## Now you are ready to try out this device

- Place the device on a flat station and sprinkle a few drops of water on the contact surface
- Turn on the power, then press the **sonic (ultrasound)** button, select low or high level
- You know ultrasound waves are generated when the water drops start to vibrate

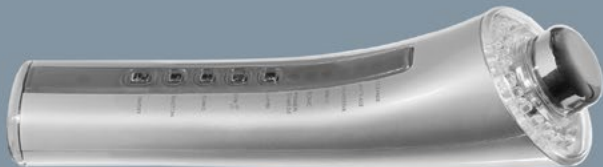

We recommend the use of our Vegas Cosmetics gel products for this treatment process in order to facilitate and maximize transmission of the ultrasound waves.


While using the ionic function, you will hear a short beep every 5 seconds. The beep will assist you in controlling the time line for the different treatment steps (for example one beep for each impurity, or 5 beeps for a dark pigment spot) and help remind you to move from point to point. You will also hear an alarm tone after 10 minutes usage time and the device will shut itself off automatically to end the treatment process. If you would like to do another treatment session, just press the on button again and choose the options you would like to use. Should you at any time during the treatment process want to stop, press the off button.





## Treatment Instructions

Apply the Vegas Cosmetics gel product onto your face (spreading it evenly without massaging it into your skin) as well as onto the contact surface. Touch the device to your face gently but ensure tight contact. If the device does not glide smoothly over your skin you may want to add some more gel (generally you will use 4 to 5 grams/0.14 to 0.17oz per application). It is very important that you always apply sufficient quantities of gel to achieve smooth movements of the ultrasound device. According to your skin type you may want to add gel several times during the treatment process.

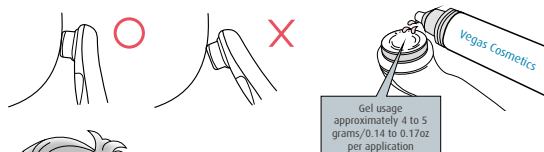


Move the device slowly within 5 seconds from your chin towards the ear. Always move upward and symmetrically on both sides of your face.

## General Care Instructions

The ultrasound device may be used up to three times per week but never longer than 30 minutes per treatment.


### First step: Deep skin cleansing process

1. Remove makeup and cleanse face thoroughly (we recommend Vegas Cosmetics Q10 Cleansing Milk or Foam Cleanser from our Vegas Products.).
2. Apply sufficient amounts of cleansing lotion or foam onto your face and the contact surface.
3. Press the **On/Off button** to turn the device on, select the **CLEANSE** function using the **COMBI** button.
4. Begin your facial massage according to the instructions in the nearby diagram regarding movements and treatment steps. Use accordingly on other treatment areas.
5. Rinse your face with water and apply a toner, i.e. Vegas Cosmetics toner. Clean device.

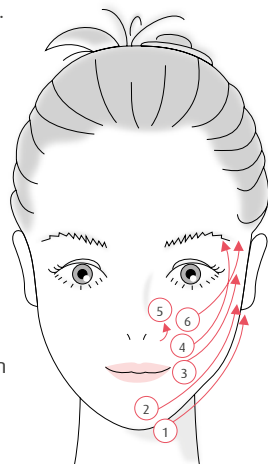

**We recommend the following Vegas Cosmetics products:** Aloe Vera foam and tonic, Q10 cleansing milk, Hyaluron Care foam and Sensitive Concept cleansing gel.





## Second step: Care for your skin and sculpt your face at the same time

1. Apply sufficient amounts of Vegas Cosmetics Hyaluronan Gel onto face and contact surface.
2. Press the **On/Off button** to turn the device on, select the **Anti-Age** function using the **COMBI** button.
3. Massage your face as described on the next page according to the diagram and move the device slowly within 5 seconds from your chin towards the ear. Use accordingly on other treatment areas.
4. Apply a moisturizer after the treatment process is completed (we recommend Vegas Cosmetics Q10 Day Cream or Night Cream, Hyaluron Day and Night Cream).

**We recommend the following Vegas Cosmetics products:** Hyaluronan Gel during the treatment process; Q10 Day Cream, Q10 Night Cream, Aloe Vera Day Cream, Aloe Vera Night Cream, or Hyaluron Day and Night Cream after the treatment process is completed.

## Third step: Photon Treatment ± LED Light Skin Care

1. Turn power on, select **photon** function and color. Red: smoothes existing wrinkles away and stimulates your body's collagen and protein production. Blue: reduces blemishes on problem skin. Green: comfortable, relaxes facial muscles.
2. Select **photon** button. Massage your skin with ultrasound and light waves, keeping the device in close contact to your face. Move the device from point to point on the areas you wish to treat.

Each area should be treated between 10 to 30 minutes.

**Hot Tip:** You may want to combine the photon treatment with the treatment process indicated in the "Second Step" by adding the desired LED light function.

## Anti Aging Care Concept

(Treating wrinkles around your eyes, mouth, nose and on your forehead)


### First step:

Cleansing your skin (same as under General Care Instructions)


### Second step:

1. Apply sufficient amounts of Vegas Cosmetics Hyaluronan Gel onto your face as well as the contact surface.
2. Press the **On/Off button** to turn the device on, select the **ANTI-AGE** function using the **COMBI** button adding Photon (red light).
3. Massage your face according to the diagram on the next page, moving the device slowly within 5 seconds from your chin towards the ear. Treat other areas accordingly.
4. After the treatment process is completed, apply a moisturizing cream (we recommend Vegas Cosmetics Hyaluron Day and Night Cream, Q10 Day Cream or Q10 Night Cream).


**We recommend the following Vegas Cosmetics products:** Vegas Cosmetics Hyaluronan Gel during the treatment process; Hyaluron Day and Night Cream or Q10 Day Cream or Q10 Night Cream for moisturizing follow-up care after the treatment process.

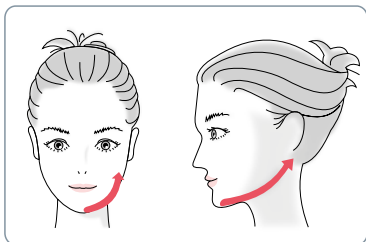


1. Always move upward  
and symmetrically

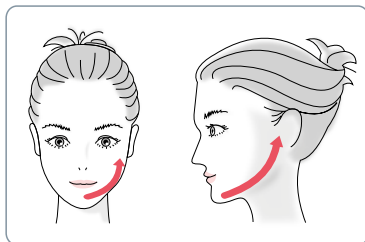


2. Use smooth, slow movements  
from chin towards your ear

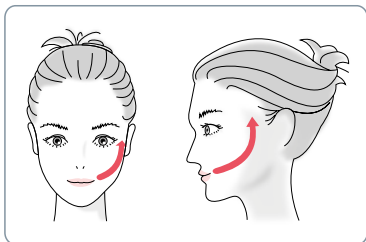


3. From the corner of your  
mouth towards the ear

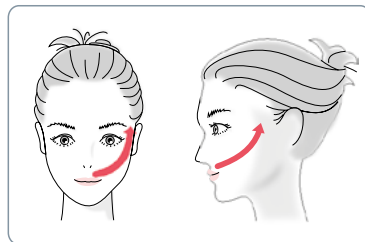

4. From your nose towards the  
end of your eyebrow

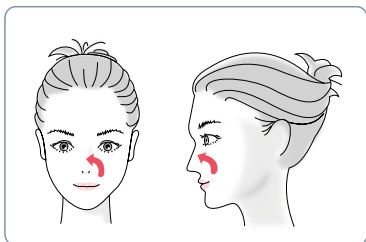

5. From under your nose to the  
bridge of your nose

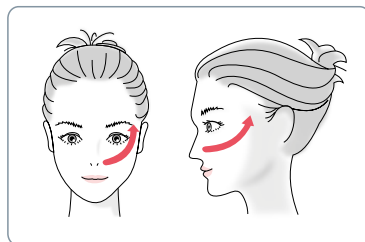

6. From your nose towards  
your temple

## Even skin tone Care Concept

(Eliminating dark pigment spots)

### First step:

Cleansing your skin (same as under General Care Instructions)

### Second step :

1. Apply Vegas Cosmetics Hyaluronan Gel onto face and contact surface. As noted previously, it is essential that you apply sufficient amounts of gel onto your skin as well as the device in order to ensure good contact as well as smooth movement.
2. Press the **On/Off button** to turn the device on, select the **Anti-Age** function using the **COMBI** button adding **Photon** (blue light).
3. Leave the device on each spot for approximately 20 to 30 seconds (wait for 5 or 6 beeps) and then move on to the next one. By keeping with this time line you will treat all spots in the appropriate manner.
4. Apply a moisturizer for post-treatment care (i.e. Vegas Cosmetics Hyaluron Day and Night Cream, Q10 Day Cream or Q10 Night Cream).

**We recommend the following Vegas Cosmetics products:** Vegas Cosmetics Hyaluronan Gel during the treatment; for moisturizing post-treatment care we recommend Hyaluron Day and Night Cream, Q10 Day Cream or Q10 Night Cream.

## Acne Care Concept

(At night)

### First step:

Cleansing your skin (same as under General Care Instructions)

### Second step:

1. Apply sufficient amounts of Vegas Cosmetics Sensitive Concept Serum onto face and contact surface.
2. Press the **On/Off button** to turn the device on, select the Sensitive function using the **COMBI** button adding **Photon** (blue light).
3. Gently move device from blemish to blemish, treating each blemish for 5 seconds (until you hear the beep) and then moving on to the next one. Proceeding in this fashion will ensure all impurities will be treated along the correct time line and, in addition, the blue LED light will help relieve infections and skin irritation.
4. Apply moisturizing cream for post-treatment care (i.e. Vegas Cosmetics 24H Face Cream).



**We recommend the following Vegas Cosmetics products:** Vegas Cosmetics Hyaluron Day & Night Cream for the initial cleansing process; Sensitive Concept Serum during the treatment process; 24H Face Cream for post-treatment care.

**Hot Tip:** If you wish, you may extend the skin care treatment by exclusively selecting the **photon** function (LED light) by itself. After the above treatment process, you may continue with the treatment process without any of the other Skin Cell Tronic functions for another 10 to 20 minutes with the blue LED light only. This will help eliminate skin irritations and infections.



### Third step: Lymphatic Detoxification

Press the **On/Off button** to turn the device on and select **SONIC (Ultrasound)**. Follow along the diagram and efficiently detoxify your lymph glands.



### Sensitive Skin Care Concept

1. Proceed according to the instructions under General Care Instructions and select **sonic (ultrasound)** function only without addition of ionic+ or ionic-.
2. You may want to add the **photon** function, selecting the red, blue or green LED light according to your choice.
3. After the treatment is finished, apply a regeneration mask (we recommend Vegas Cosmetics Aloe Vera Moisturizing Mask).
4. You may now follow your regular beauty regimen.



Note: In addition to the functions COMBI it is possible to select each function individually or combine them individually.



## Warranty Card

Product Name:

Skin Cell Tronic (Prod. code: 460)

Purchase Date:

Purchase Location:

Vegas Cosmetics Consultant:

1. In the unlikely event your device is not working properly during the first two (2) years after the purchase date, you have the peace of mind that comes with a two year warranty period.
2. Please return the device to us, together with the original sales slip, and we will repair the device at no cost to you.
3. In case the device was not used according to the manufacturer's instructions in this manual, the warranty will be void!
4. The manufacturer's warranty is explicitly void if the device was:
  - a. Opened or disassembled
  - b. Carelessly transported while in your possession
  - c. Damaged because it fell or was dropped
  - d. Intentionally damaged
  - e. Damaged by exposure to water or heat
5. The warranty card will not be replaced if lost, please keep in a safe place.



› Photon ultrasonic anti aging device

VEGAS COSMETICS GmbH  
OTTO-HAHN-STR.13  
64823 GROSS-UMSTADT

TEL.: +49 (06078) 96 81 10 | FAX: +49 (06078) 96 81 11  
[www.vegas-cosmetics.com](http://www.vegas-cosmetics.com) | [info@vegascosmetics.de](mailto:info@vegascosmetics.de)





